

300°C Laboratory Ovens

The SNOL 300LFN series from TMS Europe are a range of precision electric laboratory ovens. They are designed for low temperature thermal treatment such as drying, heating and thermal testing in an air-flow assisted environment.

Forced air convection is used to provide more effective drying and quicker heating, as well as improved temperature uniformity throughout the stainless steel lined chamber. With a temperature range of 50 to 300°C, the precision temperature control system provides good stability and uniformity for high quality results.

A digital PID temperature controller, displaying the current temperature and set-point is fitted as standard, with the option of more advanced controllers with additional functionality also available.

An independent over-temperature protection device is fitted as standard that prevents the oven exceeding its maximum safe temperature. As well as offering peace of mind, this can help meet Health & Safety and insurers' requirements, where the oven is left running unattended (e.g. overnight).

All models have a fan speed controller, as standard from TMS, allowing some adjustment of volume of air being circulated. Speed is adjustable in 10 nominal steps or the fan(s) can be turned off (only for use below 150°C). All models have a control knob to select internal or external air circulation (or a mixture as desired).

TMS Europe has these ovens manufactured with mineral insulated (MI) heating elements. These are hermetically sealed, giving longer life and more resilience to certain thermal processes, such as those involving oils, abundant moisture, carbon, etc. Subsequently the range has an increased warranty period of 2 years. Based on our many years of experience calibrating, repairing and upgrading ovens and furnaces, we offer our SNOL ovens uniquely customised with accessibility for maintenance and calibration in mind. Our sustainable design means all parts are accessible, serviceable and replaceable should the need arise. TMS have offered long-term product support on SNOL ovens for over 15 years, stocking a wide range of spares and accessories at our factory in the UK.



Our ovens are available with optional UKAS (ISO 17025) calibration of the temperature controller and sensor as a system, and we also offer multi-point thermal surveys of the chamber volume. TMS Europe Ltd is a UKAS accredited calibration laboratory No. 0461. We are ISO 17025 accredited for calibration on site and in our laboratory, as defined in our Schedule Of Accreditation (see www.tmseurope.co.uk/soa).

TMS Europe also offers models customised to meet AMS 2750G (aerospace & automotive heat treatment specification), complete with UKAS calibration and thermal survey meeting NADCAP requirements before despatch. For customisation to this specification or others please contact us.

UK Prices Exclude VAT and Delivery and are correct at the time of writing. Specification and price subject to change without notice. All trademarks acknowledged. Appearance may vary from images shown.



- *Stainless steel chamber.*
- *Air circulation fan(s) with adjustable speed.*
- *Digital PID control :- good stability, accuracy and thermal uniformity.*
- *Independent Over-Temperature Protection.*
- *Sustainable serviceable design and long-term support & repairs from TMS Europe.*
- *Many customisations and options by TMS available on short lead-times.*
- *Available from TMS's UK stock for fast delivery.*



Standard Stocked Models

Normally available from stock or on a short lead time. Options on the following pages normally take a few days longer. Please contact us for a quote with current lead-times.

Model <i>TMS Stock Code</i>	Capacity (Litres)	Chamber Size* External Size (WxDxH)	Shelves / Positions	Power / Plug	Weight	Price
SNOL 20/300 LFN <i>FCESNOL20/300LFN</i>	20 L	240 x 280 x 340mm 490 x 700 x 680mm	2 / 5	1 kW UK 13A	36 Kg	£ 1,540
SNOL 60/300 LFN <i>FCESNOL60/300LFN</i>	60 L	380 x 380 x 420mm 630 x 740 x 760mm	3 / 6	2 kW UK 13A	50 Kg	£ 1,980
SNOL 120/300 LFN <i>FCESNOL120/300LFN</i>	120 L	550 x 400 x 580mm 810 x 780 x 920mm	3 / 6	2 kW UK 13A	70 Kg	£ 2,380
SNOL 220/300 LFN <i>FCESNOL220/300LFN</i>	220 L	730 x 500 x 620mm 975 x 920 x 955mm	3 / 6	4 kW Blue 32A	102 Kg	£ 3,980

Supply Voltage: 230Vac 50Hz (single-phase). **Connector:** UK 13A = BS1363. Blue 32A = IEC60309.

Standard Non-stocked Models

Normally available on a 4-6 week lead-time. This model has double doors, hinged at each side.

Model	Capacity (Litres)	Chamber Size* External Size (WxDxH)	Shelves / Positions	Power	Weight	Price
SNOL 420/300 LFN	420 L	1000 x 500 x 860mm 1300 x 1030 x 1300mm	3 / 6	6.2kW	178kg	£ 4,880



Supply Voltage: 400Vac 50Hz (3-phase), 3P+N+E (Star with Neutral)

Connector: IEC60309 Red 5-pin 16A

* As a general rule, a gap of a minimum of 10% of each chamber dimension should be left unused on each side of the load to ensure good air flow. WxDxH = Width (Left-right) x Depth (Front-back) x Height (Top-bottom).

		SNOL 20/300	SNOL 60/300	SNOL 120/300	SNOL 220/300	SNOL 420/300
Heating Time (Minutes. Approx.)	50°C	4	5	7	3	-
	100°C	8	11	14	7	-
	200°C	19	20	27	15	-
	300°C	34	34	45	30	36
Cooling Rate (°C per Minutes. Approx.)	>250°C		3.4			
	>200°C	Contact us	2.7	Contact us	Contact us	Contact us
	>150°C		1.7			
	>100°C		1.0			
Air Changes Per Hour (With fan at full speed and vent fully open. Approx.)	50°C	4	3	2	2	-
	100°C	9	8	3,5	4	-
	200°C	14	12	6	5	-
	300°C	21	16	9	6	-
Energy Required to Maintain Temperature (KW. Approx.)	50°C	0.10	0.12	0.17	0.30	-
	100°C	0.21	0.24	0.30	0.73	-
	200°C	0.32	0.47	0.64	1.20	-
	300°C	0.49	0.78	1.01	1.60	-
Shelf Loading Weight Limit; Single Shelf / Oven Total (Kg. Non-concentrated.)		10 / 20	10 / 30	15 / 40	15 / 50	15 / 50
Packed Weight (Kg., Approx.)		56	70	100	138	228


UK Prices Exclude VAT and Delivery and are correct at the time of writing. Specification and price subject to change without notice. All trademarks acknowledged. Appearance may vary from images shown.

Temperature Controller Options

1/16 DIN size (~48x48mm) digital PID temperature controllers, with Run/Stop (Auto/Off) modes and settable heating ramp rate. With 0.1°C display resolution and featuring Autotune which can be used to optimise the control terms for the load, but which is not necessary for most applications. Optional Programmer models allow advanced timed programs (profiles) to be configured. Other instruments can also be fitted to order, providing additional options, such as audible alarms, remote communications (RS-485 etc) and data logging / recording. Please contact us for further details. The Omron E5CC and Eurotherm EPC3016 operate a buzzer that sounds briefly (user-defined) when the program finishes.

		<p>Eurotherm 3216 Programmer, 1 Program, 8 Segments</p> <p>Optional at extra cost. +£430 (Segments as Ramp+Dwell pairs)</p>
<p>Omron E5CC (Fitted as standard) Basic Timed Program: Ramp, Dwell, Stop.</p>	<p>Eurotherm 3216 Optional at extra cost. +£320 Basic Timed Program: Ramp, Dwell, Stop.</p>	<p>Eurotherm 3216 Programmer, 5 Program, 8 Segments each Optional at extra cost. +£540 (Segments as Ramp+Dwell pairs)</p>
		<p>Eurotherm EPC3016 Programmer, with RJ45 Ethernet PC Network Port,</p> <p>Optional (when OTP is also fitted) at extra cost: 1 Program, 8 Segments +£590 10 Program, 24 Segments each +£690</p> <p>Connection to a laptop or PC network via Ethernet RJ45 port on front panel. Includes PC software for easier creation, editing and backup of timed programs.</p>
<p>Eurotherm EPC3016 Programmer, 1 Program, 8 Segments</p> <p>Optional at extra cost. +£430</p> <p>Eurotherm EPC3016 Programmer, 10 Program, 24 Segments each</p> <p>Optional at extra cost. +£540</p>	<p>Eurotherm nanodac Recorder, 1/4 DIN, with RJ45 Ethernet PC Network Port, USB data transfer, Programmer option</p> <p>Contact us for full details and pricing.</p>	<p>Watlow PM Plus Programmer, 4 Programs, 10 Segments each, with Bluetooth for wireless phone app connection</p> <p>Optional (when OTP is also fitted) at extra cost. +£590 Contact us for full details.</p>
		<p>Watlow PM Plus Programmer, 4 Programs, 10 Segments each, with Bluetooth for wireless phone app connection</p> <p>Optional (when OTP is also fitted) at extra cost. +£590 Contact us for full details.</p>

UK Prices Exclude VAT and Delivery and are correct at the time of writing. Specification and price subject to change without notice. All trademarks acknowledged. Appearance may vary from images shown.

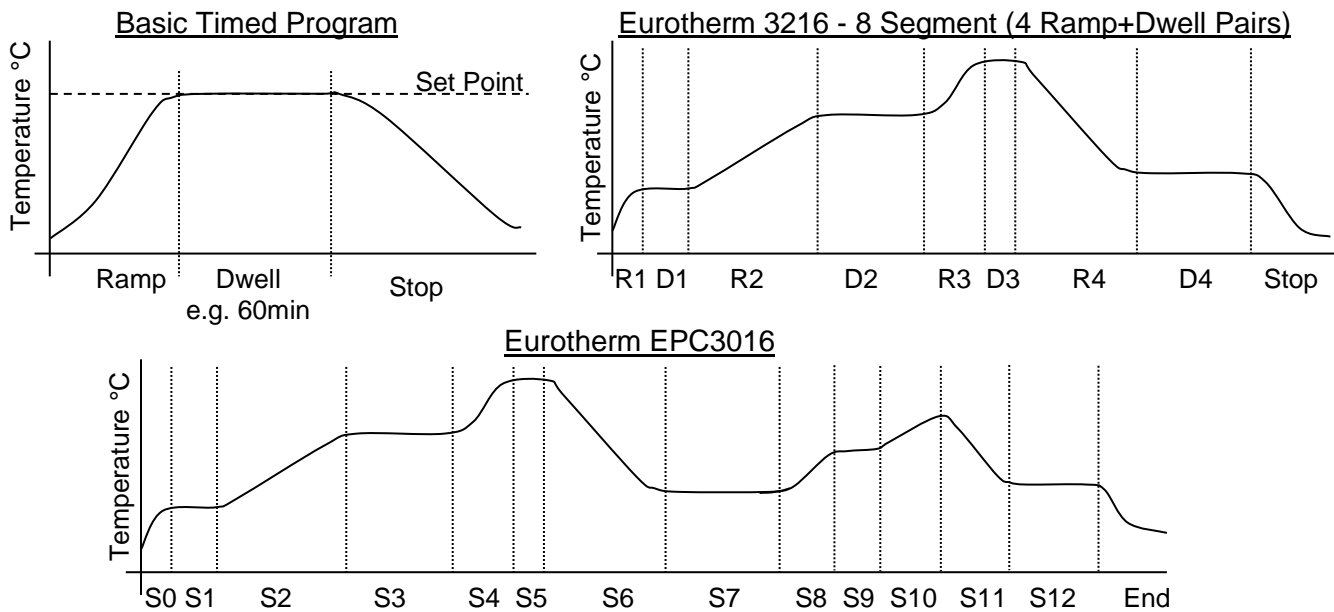

Over-temperature Protection Options

The ovens are fitted with basic over-temperature protection as standard. A user adjustable option is also offered by TMS:

<p>Basic OTP Included</p>	<p>An independent over-temperature protection thermostat device is fitted as standard that prevents the oven exceeding its maximum temperature. If tripped it, must be reset by pressing its button at the rear of the oven. As this limit is above 300°C, it won't necessary protect the user's load. <i>(Defined as DIN 12880 Class 1 category.)</i></p>
<p>OTP 2 +£270</p>	<p>Additional digital temperature limit controller mounted in the front panel. Displays its temperature reading and can be set by the user to protect their load from exceeding their desired temperature. If the oven temperature exceeds the temperature set on the over-temperature protection controller, the oven will be prevented from heating until the user resets the over-temperature protection controller, by pressing its button. <i>(Defined as DIN 12880 Class 2 category.)</i></p>

7 Day Timer +£290

A timer can be fitted to the turn on the controller at a specific time. The controller will then run at its last setting. Featuring: Daily and weekly program, Manual override, Fully automatic daylight saving time.



Accessible Thermocouple Connections

We can fit accessible thermocouple connections (miniature plug and socket) to save time with regular calibrations. If one of the following is ordered as well as OTP2 option, then connections for the OTP2 will be also be provided.

- Connections at the rear: **+£60**
- Connection in the front panel: **+£90**

Consisting of miniature size thermocouple connectors. For use by a suitably qualified person, allowing direct electrical injection onto the instrument(s) sensor input. Units with OTP 2 will have connections for both systems when this option is ordered.



UK Prices Exclude VAT and Delivery and are correct at the time of writing. Specification and price subject to change without notice. All trademarks acknowledged. Appearance may vary from images shown.

Cable Entry Ports

Option	Details	Price
22mm Survey Port (Rear)	Flange outside. In bottom right corner by default.	+£230
30mm Cable Entry Port (Rear)	Flange outside. In top left corner by default.	+£270
'50mm' Cable Entry Port (Rear)	Flange outside. In top left corner by default.	+£290
30mm Cable Entry Port (Side/Top)	Flanges outside & inside.	+£370
'50mm' Cable Entry Port (Side/Top)	Flanges outside & inside.	+£390

'50mm' sizes are nominal and in the range 47-51mm ID. Positions are as viewed from the front. Flanges protrude up to 30mm. Rear positions are limited, due to fan, vents, etc. Side positions are in the upper half of the chamber and may limit the shelf positions. Side positions can be on the left or right side, note that the door hinge is on the right. Further details are available on enquiry.

Suitable temperature rated material/bung is provided to block the Entry Port. Flanges and objects passing through the Entry Port may get hot and conduct heat outside the chamber. Excessive thermal loading or insufficient insulation in the entry port can have a negative on thermal uniformity in the chamber.



Additional Shelves

Extra Standard Shelf for SNOL 20/300 LFN	£60	Normally available from our stock.
Extra Standard Shelf for SNOL 60/300 LFN	£80	Normally available from our stock.
Extra Standard Shelf for SNOL 120/300 LFN	£100	Normally available from our stock.
Extra Standard Shelf for SNOL 220/300 LFN	£140	May be from stock, or 3-6 week lead-time.
Extra Standard Shelf for SNOL 420/300 LFN	£170	Not stocked, 3-6 week lead-time.

Calibration of SNOL Oven Systems

UKAS calibration of instrument(s) and thermocouple(s), as a system, at your choice of temperature(s) between 50 and 300°C. A calibration certificate is issued reporting the tested system's measurement at each temperature. Normally this calibration will be carried out in our laboratory prior to fitting into the oven. The typical lead-time is 4 - 10 working days. Fast turn-around may be available at additional cost.

<u>No. of Specified Temperature Points</u>	<u>Control System Only</u>	<u>Control & OTP 2 Systems</u>
1	£150	£193
2	£204	£247
3	£241	£284
4	£263	£305
5	£300	£343

A 5-point or 9-point thermal survey of the chamber volume can also be performed before despatch, please contact with your requirements for a quotation. We also provide on-site calibration services and thermal surveys.



TMS Europe Ltd is a UKAS accredited calibration laboratory No. 0461. We are ISO 17025 accredited for calibration on site and in our laboratory, as defined in our Schedule Of Accreditation (see www.tmseurope.co.uk/soa).

UK Prices Exclude VAT and Delivery and are correct at the time of writing. Specification and price subject to change without notice. All trademarks acknowledged. Appearance may vary from images shown.

Other Options

- The range is also available with stainless steel exterior (appearance may vary from unit pictured).
- Glass window (260x200mm WxH) in the door.
- Control panel at the top (for floor standing or low benches).

The lead times are typically 5 – 9 week and the price of the unit and some options are higher. Please contact us for a quote.

We may be able to assist with other bespoke requirements, please contact us with your requirements for further details.



We manufacture enhanced versions of these ovens meeting AMS 2750G and Nadcap compliant. Complete with calibration, under our ISO 17025 (UKAS) accreditation, and thermally surveyed before leaving our factory. With thermal uniformity Furnace Class 2, or better, and Instrumentation Types up to 'A'.

TMS typically holds 14 standard builds of SNOL ovens and furnaces, ready to ship or finalised within a week, in stock at our factory and warehouse (near Sheffield, UK). With additional customisations and options within 2 weeks.

Viewing and demonstration is available by arrangement.



Delivery & Shipping

These units are shipped as a palletised wooden crate, on a courier service which is normally next-day for much of the UK.

Destination	Price	Typical Transit Time
UK Mainland (England, Wales, Scottish Lowlands)	£95	1-2 days
UK Scottish Highlands, Ilse of Wight	£130	1-3 days
UK Northern Ireland & Republic of Ireland	£175	1-3 days
Other destinations (including other UK coastal islands)	Please contact us for a quote.	
Collection from our factory (Derbyshire, UK)	£30 (Weekdays, by arrangement.)	



These units are manufactured in the EU by SnolTherm (Umega Group, AB) to TMS Europe Ltd's specification. Final manufacturing, testing, localisation, customisations, addition of options and after-sale service are performed in the UK by TMS Europe Ltd. TMS was established in 1979 and has offered the SNOL range since 2008, supplying over 1,100 units to a wide range of industry sectors.



UK Prices Exclude VAT and Delivery and are correct at the time of writing. Specification and price subject to change without notice. All trademarks acknowledged. Appearance may vary from images shown.

After Sales Service & Warranty

TMS Europe has been stocking and selling the SNOL range since 2008. We hold UK stock of a wide range of spares and accessories and can offer service and repairs at our factory if required at a later date. On-site service and basic repairs may be available subject to location and the nature of the repair. We also offer a range of on-site calibration services, see our website or contact us for more details.

Each unit has a 2 year return to base warranty against manufacturing defects from the date of purchase from TMS Europe. This covers normal use of the unit in accordance with its instruction manual. It does not cover excessive 'wear and tear' or rough handling. Failures or damage resulting from the process or its by-products (e.g. leaks, spills, blockages, contamination, corrosion, etc.) cannot be considered manufacturing defects and as such are not covered by the warranty.

On receipt of the unit it is important to check for any transport damage and report it to TMS Europe and note it on the carrier's paperwork. It is recommended to keep the original wooden packaging in case the unit ever needs returning.

Under the warranty any manufacturing defects will be rectified by TMS Europe as the agent of the manufacturer at no charge. 'Return to base' means the customer is responsible for return of the unit to TMS Europe's site (Derbyshire, UK) for assessment with a view to repairing under warranty. Or, if necessary, we can provide collection at a cost, provided the unit is suitably packaged. For any work performed that is solely covered by the warranty TMS Europe will provide return shipment of the unit within the UK and Republic of Ireland at no charge. Whilst TMS Europe stocks a range of spares and aims to resolve any warranty repairs quickly, typically within 3 – 8 working days, the warranty does not guarantee this nor any provision of a loan unit while the customer's unit is with us.



TMS Europe Ltd is a UKAS accredited calibration laboratory No. 0461. We are accredited for calibration on site and in our laboratory, as defined in our Schedule Of Accreditation (see www.tmseurope.co.uk/soa).

Part of a range of thermal solutions from TMS Europe



UK Prices Exclude VAT and Delivery and are correct at the time of writing. Specification and price subject to change without notice. All trademarks acknowledged. Appearance may vary from images shown.